


KEY ACCESSORIES

PLANT GROWTH CHAMBERS

*Models: 7300-22/50/75, 7301-22/50/75, 7303-22/50/75,
7304-22/50/75, 7305-22*

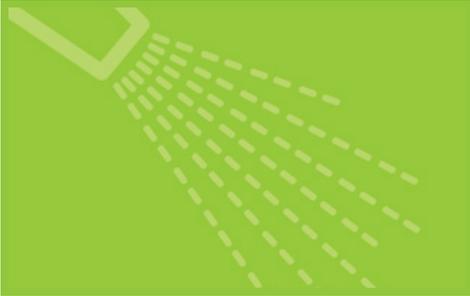

Controlled Humidity

Precise RH control without messy water pans or flooded unit bottoms. Simple atomizer-based & sensor-controlled system generates no heat and is low maintenance.



Condensate Recirculator

Provides a source of Type II equipment water for units with optional controlled RH and eliminates the need for a drain line. Includes sediment & DI filters, and a 7-gallon reservoir. Compact design fits between chambers and rolls out for easy access.



Chamber Viewing Window

Full-length sliding viewing window makes all plants visible without opening the door and disturbing growth conditions.



Ramp & Soak Programming

Ramp & soak, 5 linkable profiles, up to 50 steps total. Auto or manual ramp duration, up to 255 process loops. Easy programming, pop-up status screen included.