

KEY ACCESSORIES

PHOTOSTABILITY CHAMBERS

Models: 7540, 7545



Extended Temperature Range

Extends unit operation down to 5°C, for Q1B refrigerated product studies.




Cart

Provides a location for benchtop units, when existing counters or benches are undersized or already full. Heavy-duty aluminum construction with under-mount shelf and 4 casters.



Condensate Recirculator

Provides a source of Type II equipment water and eliminates the need for a drain line. Includes sediment & DI filters, and a 7-gallon reservoir. Compact design fits between chambers and rolls out for easy access.

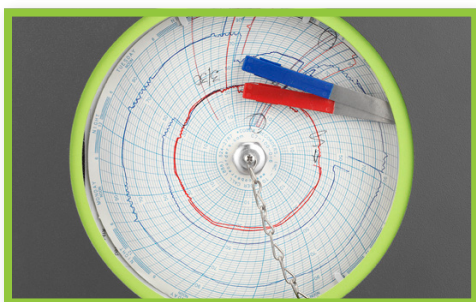

Outputs / Recorders

Circular temp/RH chart recorders in 6", 10" and 12" sizes, plus analog outputs and digital USB datalogging.